

A-dec® 200

a dec®





# The A-dec Difference.

The best way to discover what dentists need in a chair and delivery system is to ask. We did. And the A-dec 200 design reflects the way dental teams work, and purposefully supports efficiency in your practice. Everything about the A-dec 200 reflects the skill and quality of your care.



Backed by a 5-year warranty.



Designed for Durability.





**Pure light. Pure dependability.** “Good enough” just isn’t. A-dec designs and tests our lights for a minimum 20-year life.

Your equipment needs to work when you do. That’s why A-dec manages the manufacturing processes, so we can provide doctors with the most durable, reliable equipment. For example, we build and test each chair cylinder from raw material to the finished assembly, and back it with a 10-year warranty. And we manufacture our own high-quality delivery system tubing. Our production systems are consistently monitored to ensure that they meet A-dec quality standards. With the A-dec 200 solution, you get a strong, robust system designed to withstand years of constant use.



# Designed to Perform.

A-dec 200 is a system that's not only made to last, but made to grow and evolve with your practice. Detailed features such as a stain-resistant cuspidor that rotates toward the patient, a telescoping assistant's arm for close work, and a touchpad that integrates chair, light, and cuspidor functions all add up to a better workflow. For future expansion, an open platform leaves room to add or change ancillary devices. You get peak performance now, and years from now.



- Three axes light head lets you position any A-dec operatory light at just the right angle for you and your patient.
- Four preset chair positions give you added flexibility for positioning the patient.
- Ancillary accommodation allows you to add new technology as your practice evolves.
- Two-position armrests and seat articulation provide patient comfort.
- Contoured covers and minimal seams support quick and easy cleaning.

## Delivery System

Keep everything you need within reach to increase the speed of treatment.



A-dec 200



A-dec 300 Continental



A-dec 300 Traditional

## Touchpad

Fingertip access simplifies chair adjustments and provides quick access to dental light and cuspidor control. Optional deluxe touchpad adds controls for integrated clinical devices such as electric motors, cameras, and scalers.



Standard



Deluxe

## Foot Controls

Disc style with chip blower— or new A-dec lever style for precise modulation of electric and pneumatic handpieces.



Disc foot control with chip blower



Lever foot control

## Dental Lights

A-dec dental lights provide true light for accurate diagnosis, reduced eyestrain and unparalleled ergonomics.



A-dec 200 Halogen



A-dec 300 LED



A-dec 500 LED



## Monitor Mount

The monitor mount positions easily over the patient for convenient patient consultation and education.



*A-dec 200 monitor mount*



*A-dec 300 monitor mount*

## Cuspidor and Support Center

The porcelain cuspidor swivels for just the right patient position. Easy-to-fill, two-liter water bottle provides a continuous water supply through multiple procedures. The support center offers ample storage room to add ancillaries such as amalgam separators and vacuum systems.



*Cuspidor and support center*

## Assistant's Instrumentation

Three-position holder and telescoping arm conveniently positions vacuum instruments for accessibility.



*A-dec (aluminum) with touchpad*



*Durr (polymer) without touchpad*

## Formed Upholstery

A sleek vacuum-formed fit with no visible seams make cleaning and maintenance easy. See your authorized A-dec dealer for accurate color samples and the most current product information.




*sky blue*



*indigo*




*sapphire*

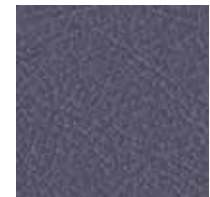


*amethyst*



*lemongrass*

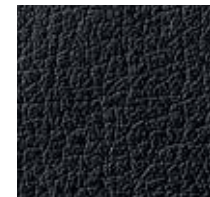

*black*

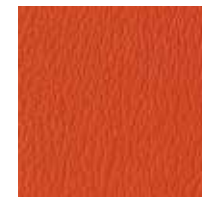

*paprika*



*tawny*

# Specifications

## A-dec 200 Dental Chair

Range of lift height	15.75" to 28.75" (400 mm to 730 mm)
Flexible backrest	2.25" (55 mm), formed
Drive system	Hydraulic base and tilt
Chair back and base limit switch	Yes, stops downward movement when obstruction is detected
Headrest	Dual-articulating (knob style)
Armrest	Pivoting
Integrated 300-watt power supply	Yes
Foot switch	Optional
Mounting capability	Floating - 2" post mount
Upholstery options	Formed

## Delivery Systems

	<b>A-dec 200</b>	<b>A-dec 300 Continental</b>	<b>A-dec 300 Traditional</b>
Control block	Standard 4-position	Standard 3-position, optional 4-position	Standard 3-position, optional 4-position
Left/right conversion	Yes	Yes	Yes
Tubing options	Vinyl or silicone	Vinyl or silicone	Vinyl or silicone
Handpiece holder positions	4	4	4
Brake handle(s)	N/A	1 standard, 2nd optional	1 standard, 2nd optional
Touchpad options	None, standard or deluxe	None, standard or deluxe	None, standard or deluxe
Quad-voltage intraoral light source	Optional	Optional	Optional
Foot control	Disc with chip blower or lever	Disc, disc with chip blower or lever	Disc, disc with chip blower or lever
Arm system	Balanced flexarm	Balanced flexarm with air brake A-dec 335 - support mount	Balanced flexarm with air brake A-dec 334 - support mount





# Building Relationships.

A-dec is headquartered in Newberg, Oregon USA, but our customers span six continents—every one except Antarctica. Schools and institutions in more than 85 countries rely on A-dec equipment to withstand the rigors of continual use day after day, year after year. Yet with A-dec, you get more than quality equipment. You get full service and support long after the sale and installation. You get local dealers and service professionals, who will respond to your needs and keep you working. And, you get the benefit of more than 50 years of experience to help your practice grow. If you ever have a question, concern or even an idea for improvement, we're just a phone call away. We're here to help you succeed, because at A-dec we don't just build equipment. We build relationships.



#### **A-dec Headquarters**

2601 Crestview Drive  
Newberg, Oregon 97132 USA  
Tel: 1.800.547.1883 within USA/CAN  
Tel: +1.503.538.7478 outside USA/CAN  
Fax: +1.503.538.0276  
a-dec.com

#### **A-dec Australia**

Tel: 1.800.225.010 within AUS  
Tel: +61.(0)2.8332.4000 outside AUS

#### **A-dec China**

Tel: 400.600.5434 within CN  
Tel: +86.571.89026088 outside CN

#### **A-dec United Kingdom**

Tel: 0800.233285 within UK  
Tel: +44.(0).24.7635.0901 outside UK



The colors and products shown may vary from the actual materials and are subject to change without notice. Please see your authorized A-dec dealer for accurate color samples and the most current product information.

All trademarks in this document are subject to the trademark rights of A-dec, its affiliates, and their collective suppliers in the United States and foreign countries.

©2021 A-dec Inc. All rights reserved.  
85.0093.00/E-only/4-21/Rev G